PARI SMARTTOUCH™ PEDIATRIC MASKS

044F5110 TODDLER LADYBUG MASK

044F5210 CHILD FROG MASK


INSTRUCTIONS FOR USE

Single patient use (one person) only.

Product is non-sterile, latex-free.

For use only with PARI Vortex® Non Electrostatic Holding Chamber.

READ ALL WARNINGS AND INSTRUCTIONS BEFORE USE.

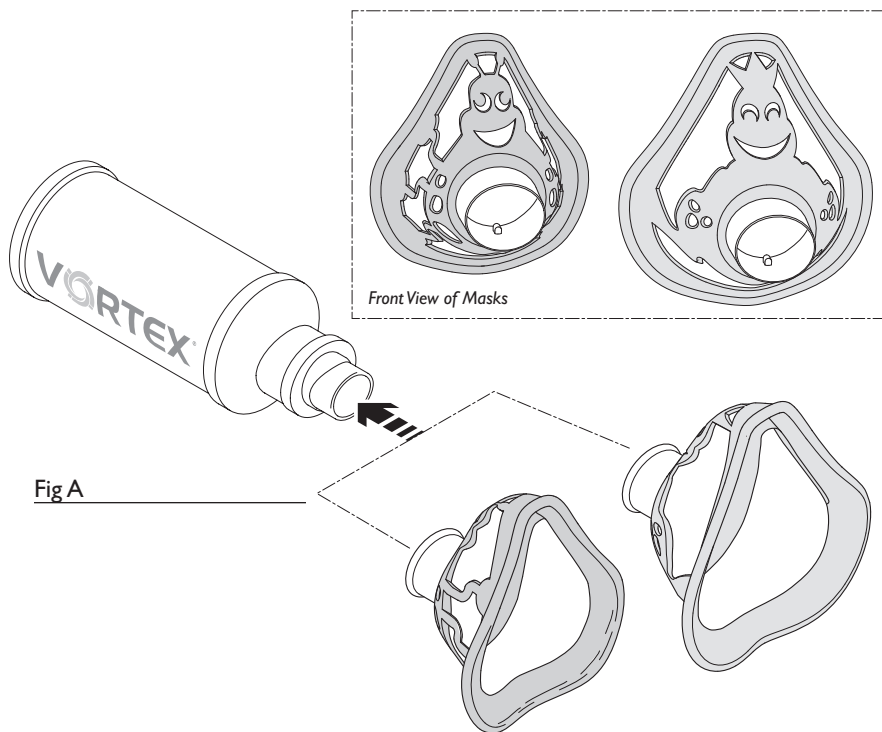Manufactured for:

PARI Respiratory Equipment, Inc.

2412 PARI Way, Midlothian, VA 23112 USA

Tel: 800.FAST.NEB (327.8632) and 804.253.7274

Fax: 800.727.4112 | Product.Info@PARI.com

Online at: www.PARI.com

Specialists in effective inhalation

© 2012 PARI Respiratory Equipment, Inc.

051D0204 Rev C 12-12

Questions? 1.800.FAST.NEB (327.8632) | PARI.com

PARI SMARTTOUCH™ PEDIATRIC MASKS

044F5110 TODDLER LADYBUG MASK

044F5210 CHILD FROG MASK

INSTRUCTIONS FOR USE

Single patient use (one person) only.

Product is non-sterile, latex-free.

For use only with PARI Vortex® Non Electrostatic Holding Chamber.

READ ALL WARNINGS AND INSTRUCTIONS BEFORE USE.

WARNING

READ ALL WARNINGS AND INSTRUCTIONS BEFORE USE. IMPROPER USE CAN RESULT IN SERIOUS OR FATAL ILLNESS OR INJURY. THESE FACE MASKS SHOULD ONLY BE USED WITH THE VORTEX® NON ELECTROSTATIC VALVED HOLDING CHAMBER. THE VORTEX® FACE MASKS SHOULD NOT BE USED WITH ANY OTHER PARI PRODUCTS SINCE THIS COULD CAUSE IMPROPER TREATMENT OR A SUFFOCATION HAZARD.

PLEASE FOLLOW STEP-BY-STEP AND REFER TO ILLUSTRATIONS ON PAGE 1.

1. **Setup:** Firmly attach your mask to the mouthpiece of your Vortex® Non Electrostatic Holding Chamber (Fig A).
2. **Taking a Treatment:** Refer to your Vortex® Non Electrostatic Holding Chamber Instructions for Use for proper Vortex® setup and operation.
3. **Cleaning:** Clean mask one of these easy ways:
 - a. Wash with warm water and liquid dish soap. Allow to dry, then store.
 - b. Place in a dishwasher basket and place on the top rack of a dishwasher. Remove and store when cycle is complete.
4. **Disinfecting:** Clean as described in number 3 above.

For low-level disinfection every other treatment day, soak mask in a solution of 1 part distilled white vinegar and 3 parts hot tap water for 1 hour. You can substitute respiratory equipment disinfectant (like Control III®) for distilled white vinegar (follow manufacturer's instructions for mixing). Rinse mask thoroughly with warm tap water and allow to dry before storing.

The mask can be disinfected by boiling in a clean pot of water for a full 10 minutes.

Hospitals and clinics can sterilize by autoclaving (a maximum of 277° F or 136° C), follow autoclave manufacturer's instructions.

Dry on a clean, lint-free cloth, then store.

Control III® is a registered trademark of Maril Products, Inc.

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Control III® es una marca comercial registrada de Maril Products, Inc.
Control III^{MD} est une marque déposée de Maril Products, Inc.